

Convenient Clinical and Technical Holter Performance.

The thin and lightweight Chroma2™ Holter recorder provides accurate pacer detection with adjustable sensitivity for the latest pacemaker technology. Add advanced analysis capabilities with easy to use, feature-rich HolterCare™ software.

FEATURES

- **Records up to 30 days of continuous, 3-channel, ECG data**
- **Data Transfer Options:** Choice of direct USB data download or removable data memory card
- **Cable Detection:** Smart cable detection between 5 and 7 lead, with display of the appropriate placement based on the cable attached
- **Direct Demographics Downloads:** Downloads patient demographics into Chroma via automatic recorder setup through direct HL7 interface
- **Adjustable Data Sample Rates:** Adjustable data sample rates from 128 to 1024 samples per second with 8, 10 and 12 bit resolution for the most demanding clinical requirements and analysis
- **Notifications:** Audible and visual patient lead disconnect and medication/activity reminders
- **Hook-up Guidance:** On-screen hook-up chart with visual guidance for proper electrode placement and lead disconnect indication
- **Multiple user-selectable languages:** Includes English, Spanish, French, Portuguese, Italian, German, Russian, Dutch, Finnish and Danish
- **Customizable:** Fully customizable recorder interface and controls
- **Convenient:** Includes neoprene pouch with belt clip

BENEFITS

- Preserves data integrity through superior ECG baseline correction for cleaner data
- Downloads ECG data in less than 15 seconds, with automated analysis in 60 seconds
- Eliminates paper diary transcriptions through easy-to-use, color coded, digital diary
- Easier to use thanks to a large, bright, color OLED display
- Confirms signal quality through visual ECG display



Chroma2: Options for Every Need

A variety of optional Chroma2 features and additional products will further enable you to maximize diagnostic and operational efficiencies.

MORE ADVANCED DATA FOR COMPLETE INSIGHT

- Multi-day Holter analysis
- Simple and reliable pacemaker analysis
- Statistical and spectral heart rate variability through HolterCare

QUICK DATA ACCESS AND TRANSFER

- EZLink allows automated transfer and analysis via the Internet

ENSURED COMPLIANCE

- HIPAA-compliant software with HL7 interface supports EHR connectivity

SINGLE PLATFORM SUITE ACCESS


CardioView Dx™ DIAGNOSTIC SUITE

- Chroma2 and HolterCare are part of CardioView™ Dx, the industry's most comprehensive single-platform, fully integrated product suite for IDTFs, hospitals and physicians
- Uniquely facilitates the management of multiple ambulatory ECG services and remote cardiac rhythm device diagnostics while reducing costs and maximizing efficiencies

The HolterCare™ Advantage

MAINTAIN COMPLETE CONTROL OVER YOUR DATA

- One mouse click provides vital information—including all arrhythmic and ischemic episodes—which can be visually confirmed and edited on a beat-by-beat basis
- HolterCare provides comprehensive arrhythmia editing techniques, including multi-level template editing, along with superimposition and page mode review
- Preview data set and set the correct channels of analysis before editing, allowing analysis of the best channels, thus yielding the best results most efficiently

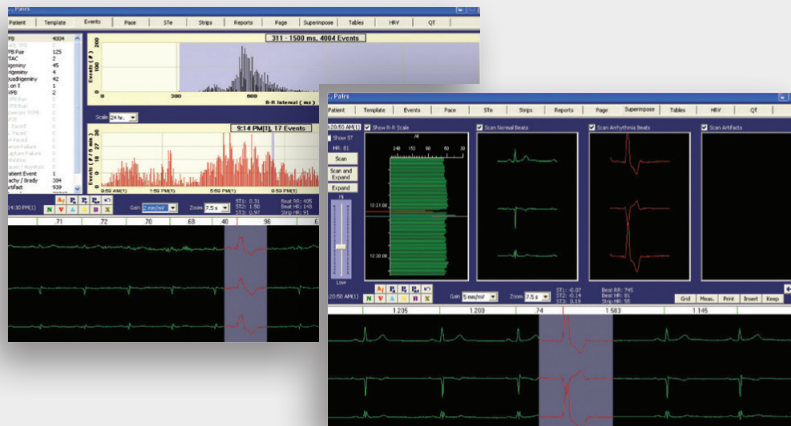

ROBUST DATA WITH MEANINGFUL PRESENTATION

- Easily locate areas of diagnostic interest through color-differentiated beats, bi-directional split-screen super imposition, page mode and colorized graphic trends
- Greatly improve efficiency with FDA XML reporting capabilities that present only the data the FDA needs to see
- Automatic atrial fibrillation identification and burden calculation utilizing industry recognized algorithm